



# High Innovation

## bms Tube Sealer TS3700 Series

Station3712

Tube Sealer-Station



Best Choice for Blood  
Component Preparation Center

Mobile3724

Tube Sealer-Mobile



Best Choice for Voluntary  
Blood Donation Center

MSS3773

Mobile Sealing Station



Best Choice for Apheresis Center

### Innovative advantages:

- Ergonomic design offers rounded edges for comfortable handling
- Sino-German R&D ensures reliable safe and efficient sealing process
- Compatible with all kinds of blood bags

Leading by Innovation

Focus • Innovation • Sustainability



# Tube Sealer 3700 series

## Station3712

### Joint R&D of Sino-German Technology team

- 8 years' joint R&D of Sino-German teams ensure fast and reliable sealing of medical PVC tube

### Modern Design

- Ergonomic design offers rounded edges for comfortable handling

### Best choice for Blood Component Preparation Center



Automatic Tube Sealer-Station3712 with sealing indicator

## Mobile3724

### Easy to carry

- Portable, battery-driven for a quick medical tube sealing
- Sufficient battery provides up to 1500 times of sealing

### Convenience to use

- Sealing handle weighs just 220g
- Up to 150cm spiral strand guarantees a flexible scope of application
- Sufficient battery life

### Best Choice for Voluntary Blood Donation Center



Automatic Tube Sealer-Mobile3724 with sealing indicator



Before sealing




During sealing



After sealing

**Usage scenarios:**

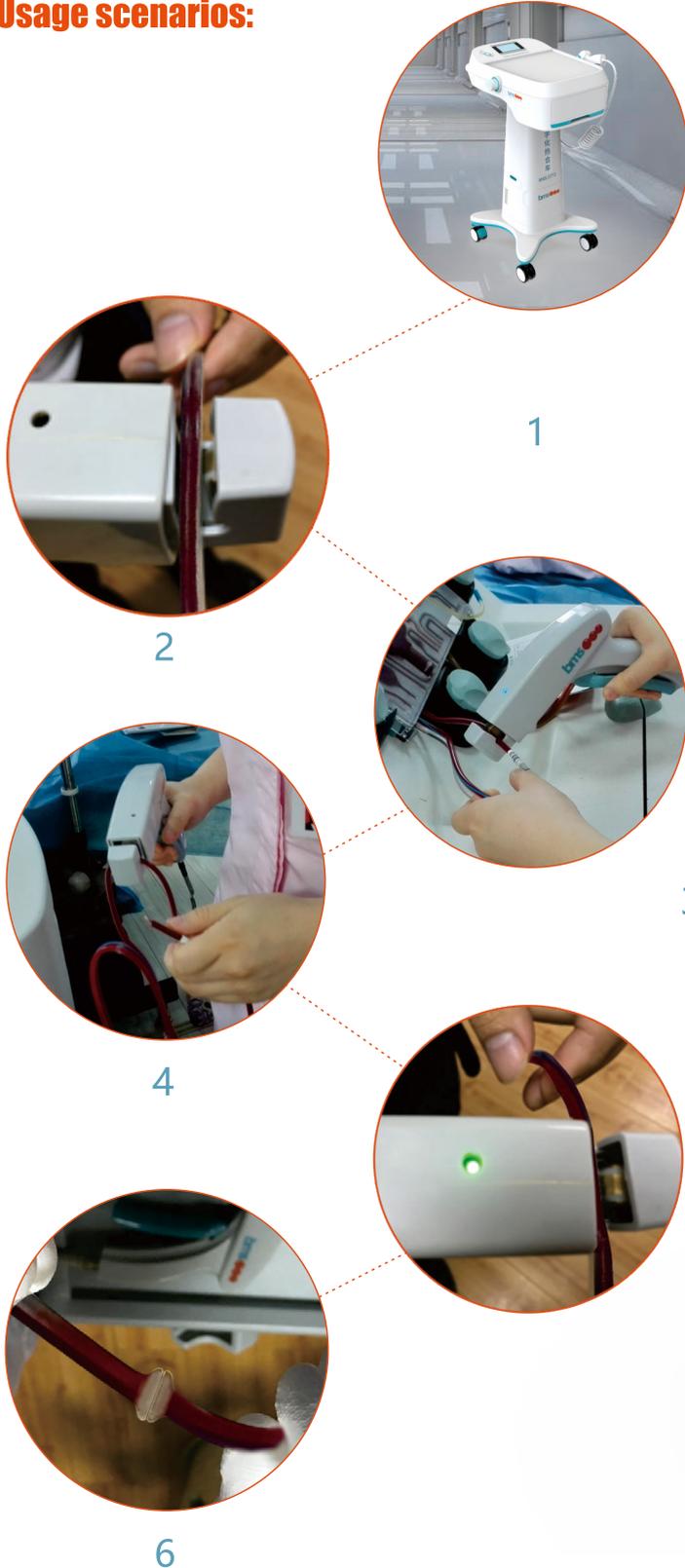

**Mobile Sealing Station MSS3773**

- Suitable for apheresis.
- Composed of Mobile Mobile3724 and a bag cart.
- Up to 20 blood bags can be carried for fixed point sealing.
- Built-in long-life lithium battery

**Best Choice for Apheresis Center**

Mobile Sealing Station MSS3773  
 with sealing indicator



## ● Technical parameters

Product Name and Model	Application	Description
Station 3712	Blood Component Preparation Center	<ol style="list-style-type: none"> <li>1.Hi frequency Dielectric heating method.</li> <li>2.Automatic sealing process to seal blood pack unit tubes.</li> <li>3.Sealing time: 0.6 to 3 sec.</li> <li>4.Capacity: upto 1000 seals per hour.</li> <li>5.Seals tubes with up to 6 mm outer diameter.</li> <li>6.Adjustable tear force.</li> </ol>
Mobile 3724	Voluntary Blood Donation Center	<ol style="list-style-type: none"> <li>1.Versatile handheld tubing sealer</li> <li>2.Hi frequency Dielectric heating method.</li> <li>3.Have a sensor automatically adjusts the sealing process based on the tubing material and diameter.</li> <li>4.Provides up to 1500 seals per charge</li> <li>5.Sealing time: 0.6 to 3 sec.</li> <li>6.Seals tubes with up to 6 mm outer diameter.</li> <li>7.LED lights for visual process control.</li> </ol>
MSS3773	Apheresis Center	<ol style="list-style-type: none"> <li>1. Integrated cart type with rollers to transfer blood bag.</li> <li>2. Multi-functional worktable reduces manual load.</li> <li>3. Built-in batteries supports 1500 times of consecutive sealing, with charging time of 4-6 hours.</li> <li>4. Indicator shows sealing state. Blue: being sealed; Green: sealing completed.</li> </ol>
Equipment Parameters		
Dimension		Station 3712 : L295MM*W180MM*H130MM Mobile 3724 : L260MM*W164MM*H110MM MSS3773 : L476MM*W369MM*H750MM Sealing gun : 180MM*145MM
Weight		Station 3712: 2.9KG Mobile 3724: 2.5KG MSS3773: 13KG
Sealing time		1 to 3 seconds

## ● Certificates and Patents

				
CE Certificate	Innovative Patent Certificates			

**bms**   
Innovation in bloodcare technology



China

Sales partner

**Wuhan bms Medicaltech Co., Ltd**

Service Hotline : +86 027 65028187

Email: info@bms-bloodcare.com

Add : 5 th Floor, C/D, B9, Hi-tech Medical Devices Park,  
818 Gaoxin Avenue, Wuhan East Lake Hi-tech Development Zone, Wuhan City, Hubei Province, China 430206