



手术床产品系列

OPERATING TABLE PRODUCT LINE


卓越品质

CH-T100 Hydraulic integrated operating table

Electric Operating Table Series



Practical electric control functions

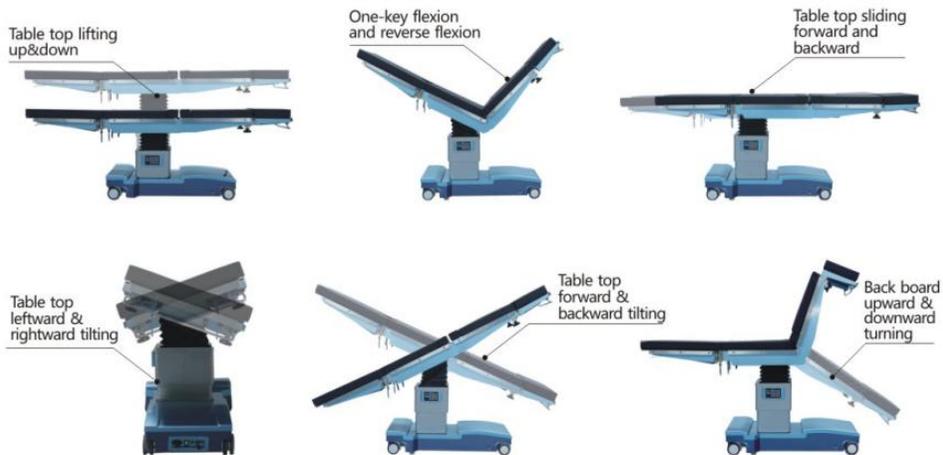


Product Features

PRODUCT CHARACTERISTICS

01. The electric hydraulic transmission structure is adopted to replace the traditional electric push rod transmission technology, making the body position more accurate, the movement speed more uniform and stable, and the performance is reliable and durable.
02. Various actions such as lifting, translation, head and foot tilting, left and right tilting, braking, etc. are driven by an independent power system.
03. The back panel adopts a carbon fiber panel, and the electric horizontal movement stroke is 320mm, which provides accurate and convenient positioning space for the C-arm, and the whole body X-ray can be performed without moving the patient.
04. The operating table can be equipped with an external waist bridge, and manual operation is convenient for more body position surgery.
05. 304 stainless steel material, stainless steel cover of the base, anti-glare treatment, reduce the eye fatigue of medical staff, and facilitate cleaning and disinfection.
06. The scientific base design leaves a wide space for the C-arm and reserves space for the feet of the medical staff to ensure the balance and stability of the operating table in achieving various actions.
07. The base of the operating table is designed with large casters (diameter $\geq 100\text{mm}$) for flexible movement. When braking, the casters rise, the bed base is firmly in contact with the ground, and has good stability.
08. The electric operating table has two modes: standard and reverse, which can realize the head and foot interchangeability.
09. There are two independent electronic operating systems, one for remote control and the other for the operating table body backup operation control panel system. The two systems operate independently to ensure that the operating table can still operate reliably when the remote control fails.
10. The operating table can be optionally equipped with high-performance rechargeable batteries of internationally renowned brands, which can meet the needs of 50-80 surgeries and ensure that the operating table can work without power. The rechargeable battery does not require maintenance and can be used for a long time. At the same time, it has AC power to provide power to ensure maximum safety.
11. The mattress is made of artificial wrapped high-density memory foam.
12. One-button flexion and anti-flexion function, using the corresponding button of

Ultra-low level hydraulic

CH-T100 electrical operating table

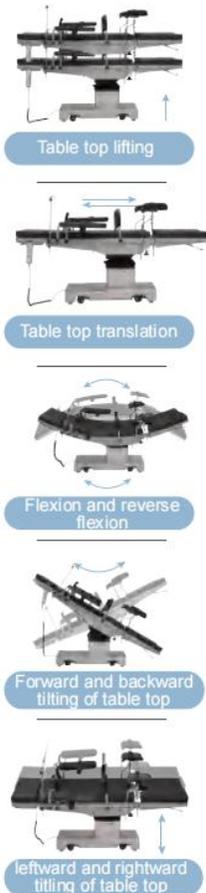

Product Features

PRODUCT CHARACTERISTICS

The electric operating table is developed by our company with the introduction of foreign advanced technology. It is a high-end comprehensive multifunctional operating table and a leading domestic product.

01. Comprehensive multifunctional operating table, suitable for orthopedics, gynecology, eye, ear, nose, throat, urology and other surgeries;
02. The main material is SUS304 stainless steel, which is sturdy and easy to use;
03. The tabletop is lifted, raised and lowered. The front, back, left, right, back and translation are electrically controlled, and the operation is smooth and convenient;
04. Low position and large lifting range, the tabletop lifting stroke is 430mm, and there is sufficient space to cooperate with the C-arm;
05. The tabletop can be translated ≥ 300 mm in all directions, and there is no blind spot for filming;
06. The electric brake is locked for easy transportation, shifting and cleaning;
07. Equipped with a built-in emergency power supply to ensure the normal use of the operating table in the event of a power outage;
08. Built-in squinting lift, cavity lifting 120mm, to meet the needs of the vertebrae, spine, gallbladder and stomach
09. One-touch flexion and anti-flexion, the bed can form a V-shaped position of 100° and an anti-V-shaped position of 220° ;
10. One-touch adjustment of automatic horizontal return is conducive to emergency treatment;
11. The bed surface uses a memory foam mattress, which is easy to clean and disinfect;
12. The manual control and column panel controller are two independent systems to control the movement of the operating table, providing smooth and precise operation. The operation is simple, safe and reliable;

Hydraulic Basis


CH-T100 electrical operating table



Table top lifting



Table top translation

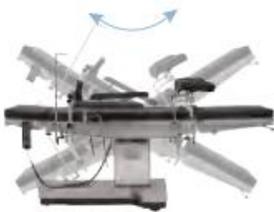

Forward and backward tilting of table top



leftward and rightward tilting of table top

Product Features

PRODUCT CHARACTERISTICS

01. The hydraulic system uses imported parts, which is stable and reliable, has a heavy load capacity and low noise, meeting the needs of different surgeries in various departments such as thoracic and abdominal surgery, neurosurgery, gynecology, urology, ENT, and orthopedics in hospitals.
02. The two independent systems of the hand controller and the column panel controller control the movement of the operating table, providing smooth and precise operation. The operation is simple, safe and reliable;
03. Modular design, the head and leg boards are detachable, and the head and feet are interchangeable. The quick-connect device greatly improves safety and accuracy, meeting the needs of surgical positions;
04. The tabletop can be translated longitudinally;
05. Professional high-strength composite bed board that can transmit X-rays is convenient for filming;
06. Special water-draining seal design, selected high-quality stainless steel material, good firmness, rust-proof, waterproof and anti-fouling;
07. Built-in waist bridge to meet the needs of gallbladder and kidney surgery;
08. One-button reset function, only one button is needed to restore the operating table to the horizontal zero position, and the bed can be restored to the horizontal position at any angle in the first time;
09. One-button flexion and anti-flexion, easily realize the waist bridge function, and meet the needs of gallbladder and kidney surgery;
10. Electric control of the operating table's rise, fall, left and right tilt, front and back tilt, translation, brake and other actions;
11. The table has a built-in waist bridge to facilitate gallbladder, kidney and other surgeries;
12. The T-shaped base design makes the operating table more stable and creates a larger activity space for medical staff;
13. Equipped with a built-in emergency power supply to ensure the normal operation of the equipment. (The built-in battery is on the base of the operating table);

High-end six functions

CH-T200 electrical operating table



Product Features

PRODUCT CHARACTERISTICS

1. The electric comprehensive operating table is used for comprehensive surgery in chest and abdominal surgery, brain surgery, ophthalmology, otolaryngology, obstetrics and gynecology, urology, orthopedics, etc., and can be used with a C-arm X-ray machine.
2. The main actions of this operating table: tabletop lifting, tabletop tilting forward and backward, tabletop tilting left and right, backboard folding, bed body translation and braking are electric devices, which are composed of a controller, electric push rods, etc.; headboard folding and legboard folding are controlled by gas springs; board lifting is mechanically controlled.
3. The tabletop can be moved longitudinally, and can be used with a C-arm for radiological examination or filming;
4. The tabletop is made of high-strength and X-ray-transmissive composite materials to meet the needs of clinical filming in hospitals.
5. The mattress is made of slow-rebound memory foam, which is comfortable and beautiful, evenly shares the load pressure, has the advantages of corrosion resistance, acid and alkali resistance, breathability, and easy cleaning. It can be naturally shaped according to the patient's body temperature to effectively prevent patients from developing bedsores.
6. The controller is equipped with a lock switch for misoperation. The handheld controller uses 24V DC voltage, which is easy to operate, safe and reliable.
7. The operating table is made of medical-grade 304 stainless steel, and the surface is specially matte treated, which is more modern and easier to clean;
8. High-strength stainless steel bottom cover has strong pressure resistance, corrosion resistance, and is easy to clean;
9. The base of the operating table is made of concave medical-grade 304 stainless steel. The base has fixed or movable functions, and the bottom has wheels, which is flexible and reliable;

Double Controlsix functions

CH-T200 electrical operating table



Product Features

PRODUCT CHARACTERISTICS

01. The overall lifting, table tilting forward and backward, table tilting left and right, backboard lifting, bed translation, and bed brake are controlled by electric devices;
02. The headboard and legboard are folded and operated by gas springs;
03. The operating table is made of medical-grade 304 stainless steel, and the surface is specially matte-treated, which is more modern and easier to clean;
04. The tabletop can be moved longitudinally, and can be used with a C-arm for X-ray examination or filming;
05. The tabletop is made of high-strength composite material that can pass X-rays (it can be upgraded to a carbon fiber bedboard), which meets the needs of clinical filming in hospitals;
06. The mattress is made of slow-rebound memory foam, which is comfortable and beautiful, evenly shares the loader's pressure, and is corrosion-resistant, acid-resistant, and alkali-resistant. Breathable and easy to clean;
07. The head and leg board are detachable, and the leg board can be manually rotated outward and folded down for easy operation, which better meets the surgical position requirements and provides medical staff with a more flexible surgical space;
08. The handheld manipulator provides smooth and precise operation. The operation is simple, safe and reliable. It can be equipped with a foot controller for more convenient operation;
09. The high-strength stainless steel bottom cover has strong pressure resistance, corrosion resistance and easy cleaning;
10. The base of the operating table has fixed or mobile functions, and the bottom is equipped with wheels, which is flexible to move and reliable to fix;
11. The table top has a built-in waist bridge to facilitate operations such as gallbladder and kidney;

Product Features

PRODUCT CHARACTERISTICS

01. The tabletop lifting, tabletop tilting, tabletop tilting, backboard folding, bed body translation, and bed body stripping are controlled by electric devices; 02. The headboard folding and legboard folding are controlled by gas springs; 03. The operating table is made of medical-grade 304 stainless steel, and the surface is specially matte-finished, which is more modern and easier to clean; 04. The tabletop can be moved longitudinally, and can be used with a C-arm for X-ray examination or filming; 05. The tabletop is made of high-strength composite materials that can pass X-rays (can be upgraded to carbon fiber bedboards), meeting the hospital's clinical filming needs; 06. The mattress is made of slow-rebound memory foam, which is comfortable and beautiful, evenly shares the load pressure, and is corrosion-resistant, acid and alkali resistant

Breathable, easy to clean;

07. The head and leg board are detachable, and the leg board can be manually rotated outward and folded down for easy operation, which better meets the surgical position requirements and provides a more flexible surgical space for medical staff;

08. The handheld manipulator provides smooth and precise operation. The operation is simple, safe and reliable. It can be equipped with a foot controller for more convenient operation;

09. The high-strength stainless steel bottom cover has strong pressure resistance, corrosion resistance and easy cleaning;

10. The base of the operating table has fixed or mobile functions, with wheels at the bottom, flexible movement and reliable fixation

11. The table has a built-in waist bridge to facilitate gallbladder, kidney and other operations;

12. The operating table is equipped with a built-in battery to ensure the normal use of the equipment in the event of a power outage.

Six functions electrical operating table



Five function operating table

CH-T200 electrical operating table



Table top lifting

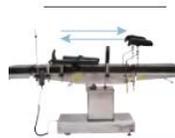


Table top translation



Forward and backward tilting of table top



leftward and rightward tilting of table top



Back lifting of table top

Product Features

PRODUCT CHARACTERISTICS

01. The tabletop is electrically operated for lifting, tilting forward and backward, tilting left and right, folding the backboard, and moving the bed;
02. The headboard and legboard are folded and operated by gas springs;
03. The lumbar board is mechanically controlled, and the lifting stroke can reach 120mm, meeting the needs of gallbladder and kidney surgery;
04. The operating table is made of medical-grade 304 stainless steel, which has strong pressure resistance, corrosion resistance, and is easy to clean;
05. The tabletop can be moved longitudinally, and can be used with a C-arm for radiographic examination or filming;
06. The tabletop is made of high-strength and X-ray-transmissive composite materials to meet the needs of clinical filming in hospitals.
07. The mattress is made of slow-rebound memory foam, which evenly shares the load pressure, is corrosion-resistant, acid-resistant, breathable, and easy to clean;
08. The head and leg board is detachable, and the separated leg board can be folded up, outward, and folded down, which is convenient for gynecological and urological operations;
09. The handheld manipulator has a smooth appearance design, conforms to ergonomics, is comfortable to hold, easy to operate, safe and reliable;
10. The high-strength stainless steel bottom cover has strong pressure resistance, corrosion resistance, and is easy to clean;
11. Mechanical brakes are convenient for fixing and moving;
12. The operating table is equipped with a built-in battery to ensure normal use of the equipment in the event of a power outage.

Product Features

PRODUCT CHARACTERISTICS

01. The tabletop is electrically operated for lifting, tilting forward and backward, tilting left and right, folding the backboard, and moving the bed;
02. The headboard and legboard are folded and operated by gas springs;
03. The lumbar board is mechanically controlled, and the lifting stroke can reach 120mm, meeting the needs of gallbladder and kidney surgery;
04. The operating table is made of medical-grade 304 stainless steel, which has strong pressure resistance, corrosion resistance, and is easy to clean;
05. The tabletop can be moved longitudinally, and can be used with a C-arm for radiographic examination or filming;
06. The tabletop is made of high-strength and X-ray-transmissive composite materials to meet the needs of clinical filming in hospitals.
07. The mattress is made of slow-rebound memory foam, which evenly shares the load pressure, is corrosion-resistant, acid-resistant, breathable, and easy to clean;
08. The head and leg board is detachable, and the separated leg board can be folded up, outward, and folded down, which is convenient for gynecological and urological operations;
09. The handheld manipulator has a smooth appearance design, conforms to ergonomics, is comfortable to hold, easy to operate, safe and reliable;
10. The high-strength stainless steel bottom cover has strong pressure resistance, corrosion resistance, and is easy to clean;
11. Mechanical brakes are convenient for fixing and moving;
12. The operating table is equipped with a built-in battery to ensure normal use of the equipment in the event of a power outage.

Five function operating table

CH-T200 electrical operating table



Table top lifting

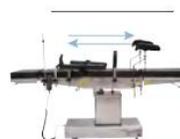


Table top translation



Forward and backward tilting of table top



leftward and rightward tilting of table top



Back lifting of table top



Four-control translation

CH-T200 electrical operating table



Table top lifting



Table top translation



Forward and backward tilting of table top



leftward and rightward tilting of table top

Product Features

PRODUCT CHARACTERISTICS

01. The tabletop lifting, tabletop tilting, tabletop tilting, and bed translation are controlled by electric devices;
02. The headboard and legboard folding are controlled by gas springs;
03. The lumbar board rises mechanically, and the lifting stroke can reach 120mm, meeting the needs of gallbladder and kidney surgery;
04. The operating table is made of medical-grade 304 stainless steel, which has strong pressure resistance, corrosion resistance, and is easy to clean;
05. The tabletop can be moved longitudinally, and can be used with a C-arm for radiographic examination or filming;
06. The tabletop is made of high-strength and X-ray-transmissive composite materials to meet the needs of clinical filming in hospitals.
07. The mattress is made of slow-rebound memory foam, which evenly shares the load pressure, is corrosion-resistant, acid-resistant, breathable, and easy to clean;
08. The head and leg board is detachable, and the separated leg board can be folded up, outward, and folded down, which is convenient for gynecological and urological operations;
09. The handheld manipulator has a smooth appearance design, conforms to ergonomics, is comfortable to hold, easy to operate, safe and reliable;
10. The high-strength stainless steel bottom cover has strong pressure resistance, corrosion resistance, and is easy to clean;
11. Mechanical brakes are convenient for fixing and moving;
12. The operating table is equipped with a built-in battery to ensure normal use of the equipment in the event of a power outage.

Four-control integration

CH-T300 electrical operating table



Table top lifting



Forward and backward tilting of table top



leftward and rightward tilting of table top



Folding of back board

Product Features

PRODUCT CHARACTERISTICS

01. Lifting, forward and backward tilting, leftward and rightward tilting of the table top, and folding of the back board are controlled by an electric device;
02. Folding of the head board and leg board is controlled by gas springs;
03. The lifting of the waist board is mechanically controlled with a lifting distance of 120mm, which meets the need of biliary and kidney surgeries;
04. The operating table is made of medical grade 304 stainless steel, which has strong pressure resistance, corrosion resistance and is easy to clean;
05. The table top can be moved longitudinally, and can be used for radiographic examination or imaging combined with a C-shaped arm;
06. The table top is made of high-strength and X-ray radiotransparent composite materials to meet the need of clinical X-ray imaging.
07. The mattress is made of slow rebound memory foam, which can evenly bear the pressure of the load and is corrosion resistant, acid and base resistant, breathable and easy to clean;
08. The head and leg boards are detachable, and the split leg boards can be folded upward, extended outward, folded downward to facilitate gynecological and urological surgeries;
09. The appearance design of the handheld controller is smooth, ergonomic, safe and reliable, and comfortable to hold;
10. The high-strength stainless-steel bottom cover has strong pressure resistance, corrosion resistance and is easy to clean;
11. The mechanical brake is convenient for movement and immobilization;
12. The operating table is equipped with a built-in battery to ensure normal use of the device in case of a power failure

Three functions

CH-T400 electrical operating table



Table top lifting



Forward and backward tilting of table top



leftward and rightward tilting of table top

Product Features

PRODUCT CHARACTERISTICS

The three-function electric comprehensive operating table designed according to the market needs using advanced technologies and meets the

needs of comprehensive surgeries in thoracic-abdominal surgery, surgery, cerebral surgery, ophthalmology, E.N.T., obstetrics and gynecology and other departments.

01. Lifting, forward and backward tilting, leftward and rightward tilting of the table top are controlled by an electric device;

02. Folding of the head board and leg board is controlled by gas springs;

03. The lifting of the waist board is mechanically controlled with a lifting distance of 120mm, which meets the need of biliary and kidney surgeries;

04. The operating table is made of medical grade 304 stainless steel, which has strong pressure resistance, corrosion resistance and is easy to clean;

05. The table top can be moved longitudinally, and can be used for radiographic examination or imaging combined with a C-shaped arm;

06. The table top is made of high-strength and X-ray radiotransparent composite materials to meet the need of clinical X-ray imaging.

07. The mattress is made of slow rebound memory foam, which can evenly bear the pressure of the load and is corrosion resistant, acid and base resistant, breathable and easy to clean;

08. The head and leg boards are detachable and the split leg boards can be folded upward, extended outward, folded downward to facilitate gynecological and urological surgeries;

09. The appearance design of the handhold controller is smooth, ergonomic, safe and reliable, and comfortable to hold;

10. The high-strength stainless-steel bottom cover has strong pressure resistance, corrosion resistance and is easy to clean;

11. The mechanical brake is convenient for movement and immobilization;

12. The operating table is equipped with a built-in battery to ensure normal use of the device in case of a power failure

Ophthalmology dedicated

CH-T400 electrical operating table



The head board is manually controlled



Electric-controlled lifting of back board



Lifting and lowering of operating table

Product Features

PRODUCT CHARACTERISTICS

01. This electric operating table is a new product developed by our company according to the clinical requirements of the hospital. It has a low-position function and is especially suitable for ophthalmology and facial plastic surgery. It is convenient for doctors to operate in a sitting position during surgery and is equipped with a hand rest to reduce the doctor's fatigue.

02. The tabletop lifting and back-lifting operation (optional) are controlled by an electric mechanism, which has the advantages of low noise, strong reliability and convenient operation.

03. The head is manually controlled, with simple operation and reasonable structure, good comprehensive mechanical properties, and scientific and reasonable design, so that the lifting mechanism can achieve stable operation, strong load-bearing capacity, low noise and other advantages.

04. The main frame is made of high-quality carbon steel, and the outer cover and lifting cover are made of 304 stainless steel, which is durable and easy to clean;

05. The mattress is made of memory foam, which effectively relieves the patient's pressure points. The material meets safety standards, the surface is antibacterial and easy to clean, and meets the disinfection requirements of the operating room;

06. The base of the operating table has fixed or mobile functions, with casters at the bottom, flexible movement and reliable fixation;

07. Equipped with accessories with various functions to meet the needs of different surgeries and expand clinical coverage;

08. The operating table is equipped with a built-in battery to ensure normal use of the equipment in the event of a power outage.



MATERNITY BED SERIES

妇产床系列

LDR electric integrated obstetric and sick table



Product Features

PRODUCT CHARACTERISTICS

An integrated table born under the guidance of a new idea of international companion delivery that integrates delivery, pre-delivery, childbirth, recuperation, gynecological surgery, gynecological examination and diagnosis.

01. Casters: 150mm TENTE centrally-controlled double-sided casters imported from Germany;
02. Hanging-type hidden dirt basin specific for LDR obstetric table (dual-dirt basin);
03. The operating system is manipulated by the control keys on the guardrails of both sides, which is convenient for parturients to adjust by themselves;
04. The host machine is a medical linear motor imported from Denmark with no need of voltage stabilizer and the operation is safe, smooth and non-noisy;
05. The adjustment including lifting and lowering of back board, knee board, the whole table and other parts are completed by a motor electric operating system and the tilting of the head board is done manually;
06. The split-type lifting guardrails make parturients safer;
07. The central-controlled brake system brings better stability for the table;
08. Starlight edged mattress or mattress with a V-shaped incision at the hip part for hospitals to choose;
09. The operating table is equipped with a built-in battery to ensure normal use of the device in case of a power failure.

Electric integrated obstetric and sick table

LDR electrical operating table



Overall lifting



Folding up of back board



Forward and backward tilting of whole table

Product Features

PRODUCT CHARACTERISTICS

The integrated maternity bed can easily realize the integration of waiting for delivery, delivery and postpartum recovery. It is suitable for delivery in various positions such as supine, sitting, side-lying, squatting, etc., providing full care, safety and reliability. 01.Steady running speed, accurate body position, reliable performance, durable:

02.ABS detachable headboard with automatic locking device, convenient for disassembly and rescue, special care and safe transportation of patients;

03.The bed body has anti-collision performance, the anti-collision corner is made of die-cast aluminum metal, equipped with anti-collision wheels, durable and can be replaced separately;

04.The guardrail is formed by polyurethane foaming mold in one time, and can be completely retracted under the bed board after being lowered. The upper part of the guardrail is in an easy-to-hold shape, which can be used as a booster for patients to stand up;

05.Equipped with accessories such as leg support and delivery handles, it is fully functional and easy to operate;

06.The central control brake silent caster is stable, safe, easy to operate, and easy to fix or move;

07.The connection between the main mattress and the foot mattress has a V-shaped incision, and the foot mattress can be folded. The mattress cover is waterproof and the filling is decompression foam, which is waterproof, soft, anti-corrosive, deodorizing and easy to clean;

08.The operating table is equipped with a built-in battery to ensure the normal use of the equipment in the event of a power outage.

Side-extending type

CH-T600 electrical operating table



Table top lifting



Upward and downward folding of back board



Forward and backward tilting of table top



Mobile auxiliary table

Product Features

PRODUCT CHARACTERISTICS

The electric side-expanding gynecological operating table is manufactured according to market demand, attracting and drawing on foreign advanced technology. The tabletop shape is designed according to human body functions, beautiful and novel; it is a multifunctional gynecological operating table for gynecological examination, diagnosis, delivery and emergency caesarean section. 01. Complete accessories and perfect functions are equipped with: backboard, armrest board, leg support, waste basin, leg board auxiliary table, anesthesia screen rack and other accessories; 02. The base is made of high-quality 304 stainless steel, which is heavy-duty, resistant to disinfection, and easy to clean, fully meeting the requirements of bed disinfection after surgery; 03. Handheld manipulator, the appearance design is smooth, ergonomic and comfortable to hold;

04. The bed is made of high-quality stainless steel, with beautiful appearance, easy to clean and disinfect;

05. The leg board auxiliary table is an axial rotation design, which can be opened 90° on the side and can be disassembled for easy storage;

06. The retractable stainless steel dirt basin can prevent the splashing of amniotic fluid during delivery and is easy to use;

07. Optional foot switch;

08. The mattress is a sponge mattress;

09. Mechanical brakes are convenient for moving and fixing the operating table;

10. The operating table is equipped with a built-in battery to ensure the normal use of the equipment in the event of a power outage.

Pull-out type

CH-T600 electrical operating table



Table top lifting



Upward and downward folding of back board



Forward and backward tilting of table top

Product Features

PRODUCT CHARACTERISTICS

The electric pull-out gynecological operating table is manufactured based on market demand, attracting and learning from foreign advanced technology. The tabletop shape is designed according to human body functions, beautiful and novel; it is a multifunctional gynecological operating table for gynecological examination, diagnosis, delivery and emergency cesarean section. 01. The tabletop is made of high-quality carbon plastic plate with excellent sealing performance. The surface is polished, impact-resistant, scratch-proof, wear-resistant, disinfection-resistant, and easy to clean, fully meeting the requirements of bed disinfection after surgery;

02. Complete accessories and perfect functions are equipped with: backboard, armrest board, leg support, dirt basin, leg board auxiliary table, anesthesia screen rack and other accessories;

03. The bed can be equipped with manual or foot control to control the movement of the operating table;

04. Built-in pull-out leg board auxiliary table, easy to store;

05. The bed body is made of high-quality 304 stainless steel, with heavy load-bearing capacity, disinfection-resistant, and easy to clean;

06. The mattress is made of memory foam mattress, which effectively relieves the patient's pressure points;

07. Equipped with mechanical brakes to facilitate the movement and fixation of the operating table;

08. The operating table is equipped with a built-in battery to ensure the normal use of the equipment in the event of a power outage.

Baby-love obstetric table

CH-T600 electrical operating table



Table top lifting



Upward and downward folding of back board



Forward and backward tilting of table top



Mobile auxiliary table



Product Features

PRODUCT CHARACTERISTICS

The electric operating table (baby-friendly delivery bed) is a multifunctional gynecological delivery bed developed and designed by our company based on market demand and the suggestions of gynecologists from various hospitals. It is a multi-purpose gynecological product suitable for gynecological examinations, diagnosis and general gynecological surgery. 01. Complete accessories and perfect functions. Equipped with: backboard, armrest board, leg support, dirt basin, leg board! Assistant table, anesthesia screen frame and other accessories;

02. Handheld manipulator, smooth appearance design, ergonomic, comfortable to hold;

03. Optional foot switch;

04. The base is made of high-quality stainless steel, with heavy load-bearing, wear-resistant and easy to clean;

05. The mattress is made of memory foam, with optional colors, antibacterial and easy to clean surface, which meets the disinfection requirements of the operating room;

06. Movable lifting auxiliary table, after the mother gives birth, it can be pushed open separately and sent to the designated area. The auxiliary table has a locking function;

07. Hidden dirt basin to prevent amniotic fluid from splashing out;

08. International standard side rails, which can be used with a variety of accessories;

09. Mechanical foot brake, convenient for fixing and moving;

10. The operating table is equipped with a built-in battery to ensure the normal use of the equipment in the event of a power outage.

Examination table

CH-T600 electrical operating table



Product Features

PRODUCT CHARACTERISTICS

This examination bed is suitable for outpatient examinations in obstetrics and gynecology, various instrument and microscope examinations, diagnosis or general surgery, etc.

01. The bed surface is streamlined and conforms to the patient's body position. Both sides of the bed are equipped with openable large arc armrests for easy adjustment of the patient's body position. The memory foam shell is soft and comfortable;

02. The handheld manipulator has a smooth appearance design, conforms to ergonomics, and is comfortable to hold;

03. The seamless high-density sponge mattress is easy to clean and corrosion-resistant, which can effectively prevent patients from developing bedsores;

04. It is equipped with a footrest with adjustable height and distance, which can be adjusted to a comfortable position according to the examination needs of the patient and the doctor;

05. A hidden dirt basin is placed under the bed, which can be easily pulled out or retracted to facilitate the treatment of sewage and debris;

06. The bed adopts advanced spraying technology, bright color, and firm attachment;

07. The bed adopts a folding structure, which is easy to clean and anti-pollution; 08. Optional plug-in auxiliary table to meet different needs;

09. Various surgical accessories to meet diversified needs;

10. The operating table is equipped with a built-in battery to ensure the normal use of the equipment in the event of a power outage.

Examination table

CH-T600 electrical operating table



Product Features

PRODUCT CHARACTERISTICS

The electric gynecological examination bed is a multi-purpose gynecological product suitable for gynecological examination, diagnosis and general surgery. This series of products has complete functions and is an ideal upgraded product for modern gynecological examination. 01. The shape of the bed is ergonomic, beautiful and comfortable, and easy to operate.

02. The lifting of the bed, the folding of the backboard up and down, and the tilting of the bed surface are all realized by electric push rods, which are flexible to operate, safe and reliable, and the angle is adjustable;

03. The height of the bed is adjustable to meet different height requirements;

04. The handheld manipulator has a smooth appearance design, is ergonomic, and comfortable to hold;

05. The mechanical foot brake is simple to operate, convenient to move and fix;

06. The retractable stainless steel dirt basin can prevent the escape of amniotic fluid during delivery and is easy to use;

07. The base is made of high-quality stainless steel, which is heavy-duty, anti-pollution, and easy to clean.

08. The table top is made of high-quality carbon plastic plate with excellent sealing performance. The surface is polished, impact-resistant, scratch-resistant, wear-resistant, disinfection-resistant, and easy to clean, fully meeting

Post-operative bed disinfection;

09. The mattress effectively relieves the patient's pressure points and prevents the occurrence of postoperative syphilis. The material meets safety standards, the surface is antibacterial and easy to clean, and meets the disinfection requirements of the operating room;

10. Configure different surgical accessories to meet various needs

Examination table

CH-T600 electrical operating table



Product Features

PRODUCT CHARACTERISTICS

1. This product is used for parturients or general obstetric surgery, especially for various instruments and mirrors for examination and diagnosis in obstetrics and gynecology and urology departments of hospitals.
2. The bed surface is divided into two parts: the head and back section and the buttocks section. It can be equipped with a plug-in auxiliary table, which is detachable and easy to operate;
3. The examination bed adopts electric adjustment and is equipped with a handheld operating system to ensure trouble-free operation. It can achieve precise linear motion and adopts a waterproof design to prevent sewage, dirt, and static electricity, so that it can be used flexibly and freely for many years;
4. It is suitable for obstetrics and gynecology outpatient examinations, various instrument mirrors, diagnosis or general surgery. The bed surface design is streamlined and conforms to the patient's position;
5. Equipped with a footrest with adjustable height and distance, it can be adjusted to a comfortable position according to the examination needs of the patient and the doctor;
6. A hidden dirt basin is placed under the bed, which can be easily pulled out or retracted to facilitate the disposal of sewage and debris;
7. High-strength spray-coated bed, easy to clean, anti-pollution, folding structure, high load-bearing and high stability;
8. Seamless high-density mattress, easy to clean, corrosion-resistant, not easy to break, antibacterial, anti-static, can effectively prevent patients from bedsores;
9. Built-in emergency power supply, can automatically jump to working state when power is off

Orthopedic traction rack



Product Features

PRODUCT CHARACTERISTICS

01. The main metal parts are made of high-quality stainless steel and are durable;
02. The pelvic frame has three optional plug-in positions: left, middle, and right. When performing hip joint fluoroscopic surgery, it can avoid the obstruction of the traction bracket and allow the doctor to perform surgery closer to the patient on both sides;
03. It is easy to use with any position of the C-arm, without obstacles, dead angles, and full fluoroscopy;
04. The pelvic frame on the lifting column can also be plugged into the side-lying frame or knee joint frame for use. The rising and falling stroke is 30 cm, which is enough to adapt to the differences in the size of bones of men, women, old and young;
05. The traction rod can change the angle greatly and has a telescopic design. In addition to the 360° rotation of the ankle fixator itself, it can also be removed to install other accessories, which is flexible in application;



Product Features

PRODUCT CHARACTERISTICS

Product performance structure and composition This product consists of a connecting fixture, a perineal column, an abduction rod, a traction device, a foot plate fixing rod, a side rail, a support column, a steering joint, a frame, a traction foot plate, and a traction foot plate direction clamp.

Product Role

Progress of operation.

During the operation, the use of a traction frame allows the patient to easily adjust the traction device under the fluoroscopy of the C-arm machine to achieve a good reduction effect, and it is not easy to change the body position, which significantly helps the operation progress of the operation.